

MEDICAL CABINET MLK-65

The MLK-65 medicine cabinet, based on the Merianto ISO module system, has been developed taking into consideration the needs of medical employees for a universal medicine cabinet. The lockable medicine cabinet has two parts and contains a separate lockable metal safe. The contents of the cabinet can be custom assembled using modular ISO system baskets of 600×400 mm or glass shelves.




The medicine cabinet is produced from quality melamine. The selection for the facade includes a standard door (left/right), as well as a space-saving roller door. The door opens to 220°.

The load capacity of the shelves, drawers and baskets is 20 kg.


There are many different ways of assembling the cabinet, as the cabinet's interior components – side panels, baskets, drawers, and shelves – can be assembled to specific requirements.

Standard options include five different finishing's:

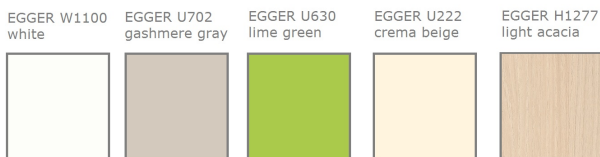

It is also possible to select from among tens of different colours and finishes from the catalogues of melamine producers, the availability of which must be verified at the time of ordering.

- The standard solution of the MLK-65 medicine cabinet includes:
- 18 mm quality melamine shell and doors. The selection of materials also includes an antibacterial surface finish.
- Pull-out auxiliary surface on ball runners, soft closing and lockable into the open position.
- There are four adjustable glass shelves in the upper part of the cabinet that can be replaced with ISO baskets if required.
- The upper part of the cabinet contains a lockable metal cabinet, white RAL9010.
- The upper part of the cabinet has sensor-controlled LED lighting.
- The interior surface of the door of the upper part has Plexiglas height adjustable pockets, with dividers.
- The lower part of the cabinet contains ISO 600×400 mm baskets, with dividers and label holders. The sets of baskets can be selected from within SET1 – SET4.
- Locks on all doors. The locks can be serialized as desired so that they are compatible with other cabinets in the same room, floor or similar area. The base selection includes a mechanical furniture lock, but an electronic code lock or chip card lock can also be selected.
- The silicon magnets on the doors ensure the secure and soft closing of the doors.
- Strong metal foot frame with adjustable foot ends, colour white RAL9010.
- Handles are 128 mm long and have an antibacterial surface finish.

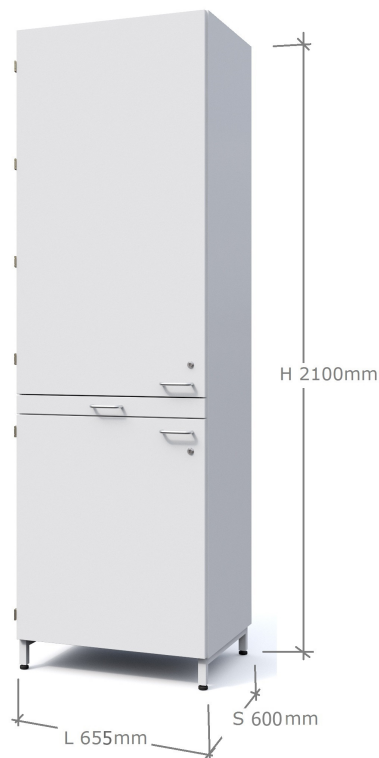

MEDICAL CABINET MLK-65

Lower cabinet basket selection (basket height mm)



SET 1

50
+
100
+
100
+
100
+
100



SET 4

100
+
100
+
100
+
200



SET 2

50
+
100
+
200
+
200



SET 5

100
+
100
+
100
+
100
+
100



SET 3

50
+
50
+
50
+
50
+
50
+
100
+
100



SET 6

200
+
200
+
200